ISNATCH

Cod. 64.2506.50

TERMOMETRO IR PER TEMPERATURA CORPOREA (IT)

EN - Infrared thermometer for body temperature and surface
FR - Thermomètre IR pour la température
DE - IR-Thermometer für Körpertemperatur oder Oberflächen
ES - Termómetro IR para temperatura del cuerpo o de superficies
PT - Termómetro IV para temperatura corporal ou superfícies

CE

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(EN) User manual - (FR) Manuel d'utilisation - (DE) Bedienungsanleitung
(ES) Manual del usuario - (PT) Manual do usuário

KON.EL.CO S.p.A. - Piazza Don E. Mapelli, 75 - 20099 - Sesto S. Giovanni - MI - Italy

MADE IN CHINA

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EN - Due to the continuous evolution of the products, the features and design of this model may change without notice.

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EN - The crossed-out wheeled bin symbol shown on the device indicates that the product at the end of its useful life must be collected separately from other waste. Disposal information is available on the website www.gbconline.it.

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EN - Device for domestic use, suitable for checking body temperature for preventive purposes. Consult a doctor to carry out any diagnosis.

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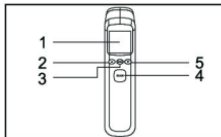
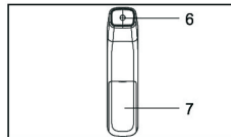
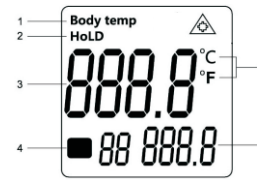
EN

SYSTEM	CONTACTLESS
Operating distance	0 ÷ 20 mm
Body measurement range	"Normal range: 32,2°C ÷ 42,2°C ±0,2°C It shows "LO" when lower than 32°C, It shows "HI" when higher than 42°C"
Operating temperature	10°C ÷ 40°C (50°F ÷ 104°F) Humidity ≤ 85%
Storage temperature	20°C ÷ 55°C (-4°F ÷ 131°F) Humidity ≤ 85%
Power supply	3V 1,5W (2x 1,5V AAA batteries)
Power Dissipation	≤ 30 mW
Battery level indication	Yes
Memory	32 groups
Display unit of measure	(°C) / (°F)
Automatic shutdown	After 40 seconds
Dimensions	166 x 39 x 40 mm
Net weight	68,6 g (batteries not included)

(IT) BATTERIE NON IN DOTAZIONE - (EN) BATTERIES NOT SUPPLIED - (FR) BATTERIES NON FOURNIES - (DE) BATTERIEN NICHT ENTHALTEN - (ES) BATERÍAS NO SUMINISTRADAS - (PT) BATERIAS NÃO FORNECIDAS

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EN - Use

Unit of measure

Operation: Press the ▼ button to change the unit of measurement and choose Celsius °C or Fahrenheit

Mode Selection

Operation: In standby mode, press ▲ button to choose "Body temp" mode for body measurements or "HoLD" for the objects. (if the display indicates "LOG" try again).

Measurement history

Operation: In standby mode press MODE for 3 seconds when LOG appears. Press ▼ or ▲ button to view the data. Press MODE again for 3 seconds to exit.

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EN - Troubleshooting guide

Indication on display	Explanation
HI	What is being measured is above the maximum reading limit of the instrument (over 42°C)
Lo	What is being measured is below the minimum reading limit of the instrument (below 32°C)
	Replace batteries

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EN - WARNING

Keep the instrument clean. Clean the instrument with a soft cloth, do not use chemicals to clean the instrument. Keep the instrument away from water (do not immerse in liquids) or heat sources including direct exposure to sunlight. Do not place batteries with inverted positive and negative polarity. Remove the batteries from the instrument if you will not use it for a long period (over 3 months). Do not use it if there are signs of damage to the sensor. Do not expose the instrument to fire as it may explode.

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